


SPECIFICATION

Infrared Thermometer / FC-IR202



- ① LCD display screen
- ② Measure button
- ③ Mem (Memory Recall)
- ④ Mute-unmute button (or °C/°F conversion)
- ⑤ Sensor
- ⑥ Rating label
- ⑦ Battery cover

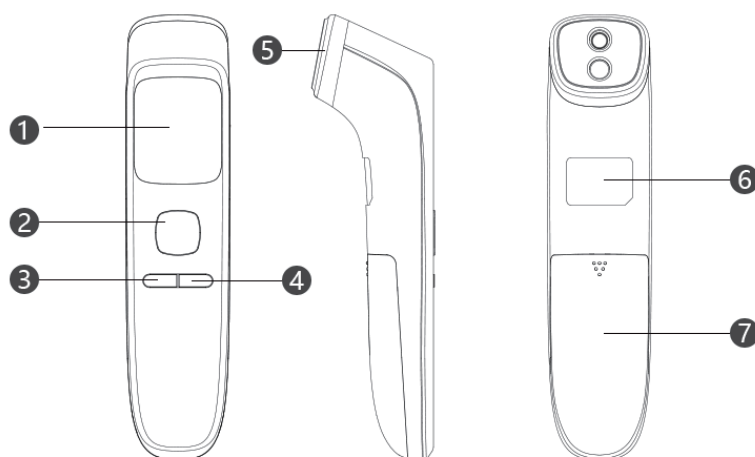

Forehead Measure



Baby mode



Fever Warning



°F/°C
Switchable



Clinically Tested



Object/Room
Measure

FC-IR202



FEATURES

- One Second Reading
- Tri-color Backlit
- Accu-sensor Technology
- 35 Groups Memory Function

● FEATURES

- Smart and beautiful.
- Baby mode switchable, with more accurate measurements for kids.
- Instant reading; Quick measurement, get reading in 1 second.
- Adopt the latest algorithm with optimised temperature calibration function.
- Accurate and reliable; $\pm 0.2^{\circ}\text{C}$ High accuracy.
- Tri-color back light, display green, orange & red with fever alarm function
- Memory function; Memory recall of 35 reading.
- Multi function design, can measure forehead/ear, room, object temperature.
- Switching between mute and un-mute.
- Switching between Fahrenheit and Celsius.
- Automatic shutdown and power saving.

● SPECIFICATIONS

| Model No. | FC-IR202 |
|---------------------------|--|
| Measurement range | Forehead: 32.0°C - 42.9°C (89.6°F - 109.2°F) Object: 0°C - 100°C (32°F - 212°F) |
| Measuring accuracy | Forehead mode: $\pm 0.2^{\circ}\text{C}$ / $\pm 0.4^{\circ}\text{F}$ Object mode: $\pm 1.0^{\circ}\text{C}$ / 1.8°F |
| Display | LCD backlight |
| Automatic shutdown | 10s \pm 1s |
| Memory | Memory 35 groups of measured temperature. |
| Indicator and Fever alarm | Normal green light: $32^{\circ}\text{C} \leq T \leq 37.3^{\circ}\text{C}$ (89.6°F - 99.2°F) Low fever orange light: $37.4^{\circ}\text{C} \leq T \leq 37.9^{\circ}\text{C}$ (99.3°F - 100.3°F) High fever red light: $38^{\circ}\text{C} \leq T \leq 42.9^{\circ}\text{C}$ (100.4°F - 109.2°F) |
| Operational conditions | Temperature: 10°C - 40°C (50°F - 104°F) Humidity: 15-95%RH, non-condensing Atmospheric pressure: 86-106 kPa |
| Battery | 2*AAA, can be used for more than 3000 times |
| Weight | 66g(without battery); 81g(with battery) |
| Product dimension | 36*42*153.5mm |
| Packing size | 46.5*46.5*169mm |
| Full set weight | 115g(including products/bags/instructions/battery) |